11-90-001



Veloped Cane holder – Assembly Instruction


2-piece cane holder for the Veloped that also works as a safety strap when the Veloped is folded for transport or storage. The cane holder fits all Veloped models.

Scan the QR-code with your smart phone or device to view video instructions.



(1) Attach the rear holder.

Wrap the rear holder around the upper part of the rear vertical frame tube, just below the height-adjustment quick-release.

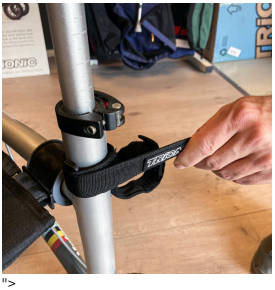

(2) Pull the strap tight.

Wrap the strap around the tube and attach it with the Velcro. Let the excessive length of the strap hang freely.

">



(3) Adjust the rear loop.

The size of the rear loop can be adjusted by loosening the Velcro closure. You may need to gently squeeze the loop to get a rounded/circular shape.



4 Attach the front holder.

Wrap the front holder around the upper frame tube, just behind the basket rail attachment. Make sure that the anti-slip patch is facing the frame tube.



5 Adjust the front loop.

The size of the front loop can be adjusted by loosening the Velcro closure. You may need to gently squeeze the loop to get a rounded/circular shape.



6 Attach the walking stick.

Attach the walking stick by pushing it from the rear, first through the rear and then through the front holder.



7 Secure the frame.

You can use the cane holder strap to secure the frame when it is folded, when transporting or storing your Veloped.